

Fall Management Technology



krug HEALTHCARE 

Fall Management Technology

Krug's proprietary fall management technology is designed to support a facility's fall management program. The built-in imperceptible sensor is intended to minimize the amount of unassisted exits and patient falls, by notifying caregivers when a patient deemed a fall risk attempts an unassisted exit from a recliner.

Since patients are often unattended while in their room, there's a much higher chance of fall occurring in this space than in other parts of the facility. Krug's simple fall management technology helps minimize the risk of falls when patient's attempt to make an unattended exits from the chair in their room.





Benefits

Falls & Injuries: The primary benefit is real-time alerts when someone tries to exit the chair, enabling staff to intervene before a fall occurs, which is crucial as falls often lead to injuries like hip fractures or head trauma.

Provides Peace of Mind: Caregivers and family members gain reassurance that their loved ones are safer, reducing anxiety about accidents happening when they're not present.

Supports Safe Mobility: They help manage risky movements, like getting in and out of the chair, by alerting staff to assist with transfers, preventing falls from happening during these transitions.

Discrete Monitoring: A modern less intrusive system offering a more dignified experience.

Reduces Wandering: For individuals with cognitive decline, these alarms can also help prevent wandering by alerting staff to attempts to leave a safe area.



Ascend Recliners



Ascend Exam Recliner



Ascend 3-Position Recliner



Ascend Sleep Recliner



Ascend Recliner Plus

Jordan Recliners & Gliders



Jordan 3-Position Recliner



Jordan Sleep Recliner



Jordan Recliner Plus



Jordan Compact Recliner



Jordan Active Patient Recliner




Jordan Glider




Integrated Wire Management



Easy Access Monitor

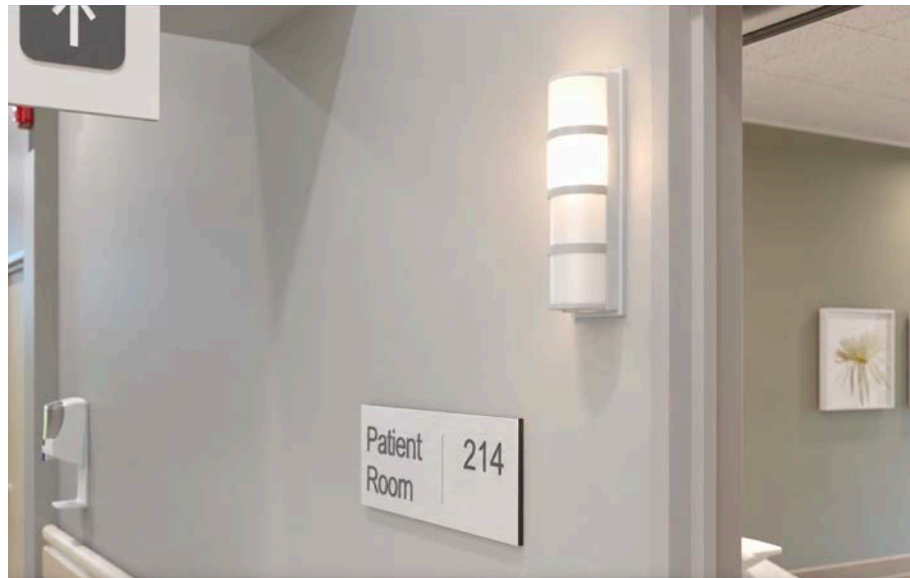

Wireless Patient Room Alarm

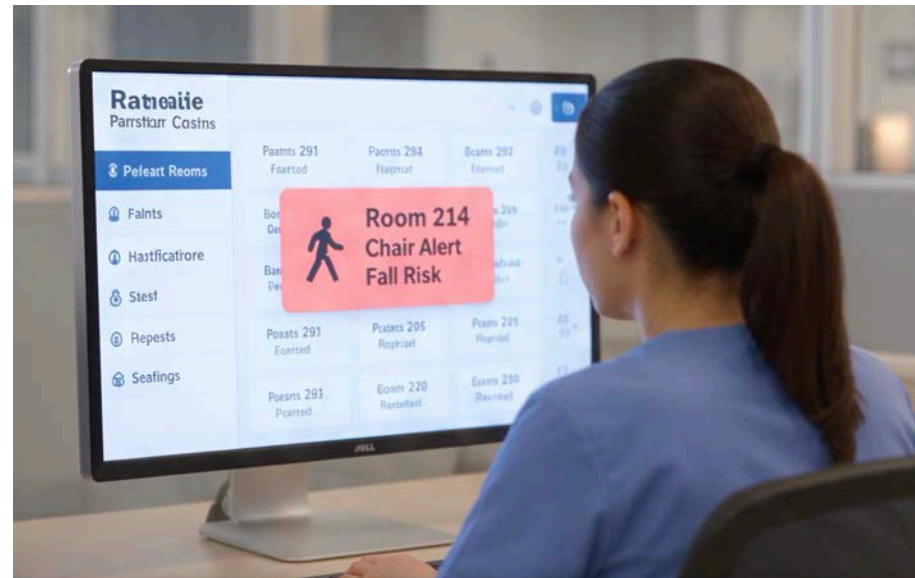

Nurse's Station Alert Software



Frequently Asked Questions

How does the Krug Fall Management sensor work?

The sensor is built into the seat of the chair. There is a cord on the sensor that attaches to the monitor on the back of the recliner. There is a second cord that connects into the monitor and the hospital's nurse call system.

What nurse call systems are compatible with Krug's Fall Management sensor?

Krug's fall management sensor is compatible with virtually all existing nurse call systems- whether that be a system that has been installed for decades to one that was installed yesterday including these systems-

What is the warranty on the Krug Fall Management sensor?

The sensor has a 2 year warranty

Is the Fall Management sensor compatible with wireless nurse call systems?

We are working on a wireless option- not yet available

How is the monitor powered?

Conventional Corded Power or Battery Power

How do you replace the sensor if it becomes damaged?

The seat will need to be replaced. Removable seat covers can be specified for ease of replacement.

How do you replace the monitor if it becomes damaged

A replacement monitor will fit into the existing holder.

Nurse Call Systems

ASCOM®*/GE

18 pin connector
Telligence™

AUSTCO®

8 pin DIN connector
MediCom
Tacera

CRITICAL ALERT®*

8 pin modular connector
SPS101
System 200, 1000/3000

DUKANE

17 pin connector
ProCare 2000®
System 2070

DUKANE

18 pin connector
ProCare 4000
ProCare 6000
System 2070

FISHER BERKELEY

15 pin connector
HP Series
MMC Series

HILLROM®

20 pin connector
COMposer®
NaviCare®

HILLROM®

37 pin connector
DB37
Headwall

INTERCALL

8 pin DIN connector
Legend II, Equinox II,
Spectrum II,
Ultra 6000, 6900, 9000

JERON®

8 pin DIN connector
Pro-Alert™ 640
Provider® 680 and 790
Spectrum® Pro IR

RAULAND

8 pin DIN connector
Responder III+, IV
Responder 4000, 5, 5000
Responder Enterprise

SIMPLEXGRINNELL®

15 pin connector
EZCare™

TEKTONE®

8 pin DIN connector
Tek-CARE® NC300 II
Tek-CARE 400

TEKTONE®

8 pin modular connector
Tek-CARE NC300 II

WEST-COM

8 pin modular connector
1000, 8000
Endeavor
Marc II

WEST-COM

10 pin modular connector
Odyssey,
Odyssey Plus, Novus

ZETTLER

9 pin connector
Sentinel 400™,
Sentinel 500™

Krug Environmental

The use of recycled and recyclable materials, combined with environmentally-friendly production processes and packaging, have made Krug's Healthcare products level™-certified, and contributors to LEED™-driven specifications. Krug Healthcare products are certified by the Green-guard™ standard for indoor air emissions, and their unique capabilities for extending product life makes them a responsible environmental choice - and one that will provide exceptional return on investment.



All products are air emissions certified and are available as FSC® certified



krug HEALTHCARE 

solutions@krug.ca | www.krug.ca | 1.888.578.KRUG